

Detection of interproximal caries *in vitro* using The Canary System[®]

B. Wong¹, S. H. Abrams¹, C. Tasevski², K. Sivagurunathan¹, R. J. Jeon¹, J. D. Silvertown¹,
W. M. P. Hellen², G. I. Elman², B. T. Amaechi³

¹Quantum Dental Technologies, Toronto, ON, Canada

²Cliffcrest Dental Office, Scarborough, ON, Canada

³Comprehensive Dentistry, University of Texas Health Science Center at San Antonio, San Antonio, TX

Introduction

- Detection of non-cavitated interproximal caries is of great importance because disease progression may be halted at this stage, remineralized or minimally restored, thereby preserving natural tooth structure.
- Visual and tactile methods of interproximal caries detection is challenging due to inaccessibility.
- Interproximal lesions at contact point can be difficult to identify on radiographs.
- Therefore, an adjunct method that can assist clinicians in the detection and quantification of interproximal caries would be of value in the prevention and management of dental caries.



Objective

- To evaluate the ability of The Canary System, International Caries Detection and Assessment System (ICDAS II), and radiographic examination to detect natural interproximal decay *in vitro*.

Caries Codes

- 0 = Sound tooth surface
- 1 = First visual change in enamel
- 2 = Distinct visual change in enamel
- 3 = Enamel breakdown, no dentine visible
- 4 = Underlying dentinal shadow (not cavitated into dentine)
- 5 = Distinct cavity with visible dentine
- 6 = Extensive distinct cavity with visible dentine



**Visual Examination
(ICDAS II)**

**Radiographic
Examination**



The Canary System®


The Science Behind The Canary System®

➤ Energy Conversion Technology

- Pulses (2 Hz) of laser light (660 nm.) are shone on the tooth surface for 5 seconds.
- Tooth glows (Luminescence, LUM) and releases heat (Photo-Thermal Radiometry, PTR).
- Canary algorithm combines detected signals to create a **Canary Number**, which reflects the tooth's state of mineralization and crystallization.
- Dental caries affect PTR-LUM signals.
- Detects 50 micron lesion up to 5 mm below the surface.



The Canary Scale

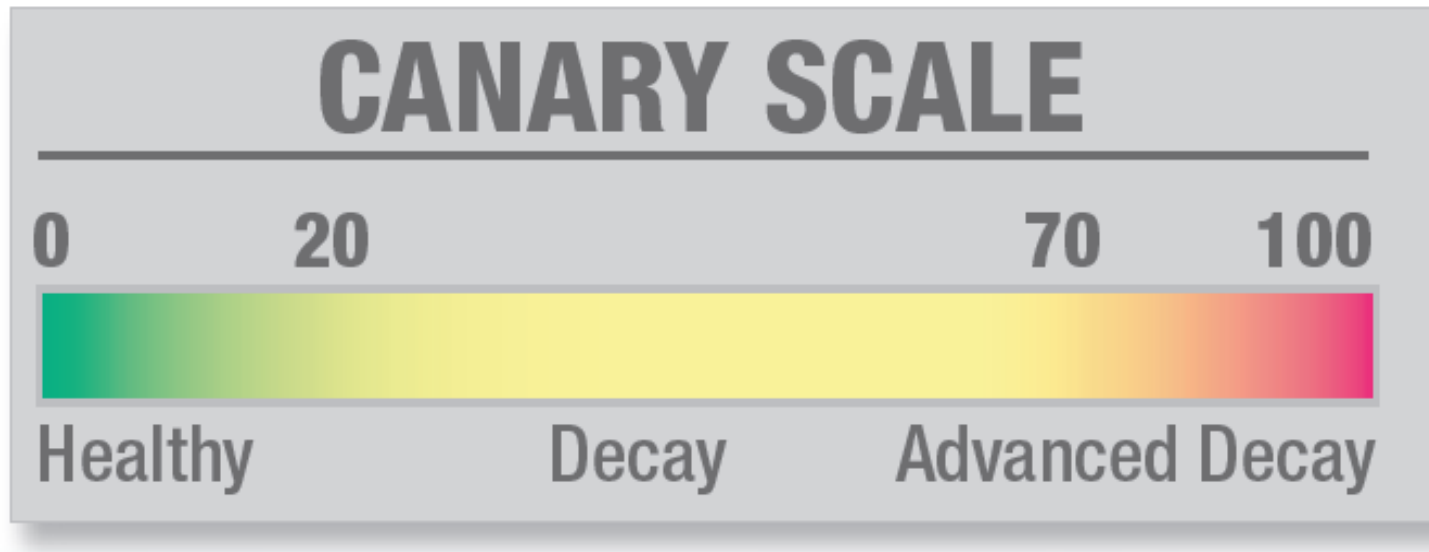


Previous Studies

- Demonstrated the ability of the core technology of The Canary System, called photothermal radiometry-luminescence (or PTR-LUM) technology, to detect very early artificially demineralized interproximal lesions.
 - Jeon RJ et al. J Bio Optics. 2007;12(3):034028.
 - Jeon RJ et al. Caries Res. 2006;40:348.
- PTR-LUM corroborated with μ -CT, TMR, SEM and PLM.
 - Jeon RJ et al. J Bio Optics. 2007;12(3):034028.
 - Mandelis A et al. Eur Phys J. Special Topics. 2008;153:467-469.
 - Jeon RJ et al. J Bio Optics. 2008;13(3);034025.
 - Matvienko A et al. Proc. SPIE BIOS. 2009;7166 (12);71660C1-12.
 - Wong et al. J Dent Res. 92 (Spec. Iss. A), 7, 2013.



Materials and Methods


- Twenty interproximal surfaces of ten pairs of extracted permanent human teeth were examined.

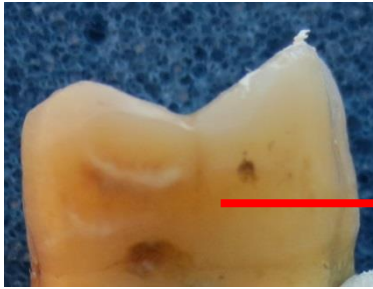

Visually carious
interproximal surface
of Tooth #1



Adapted from: Buchalla W et al.
Caries Research. 2002;36:320–326.

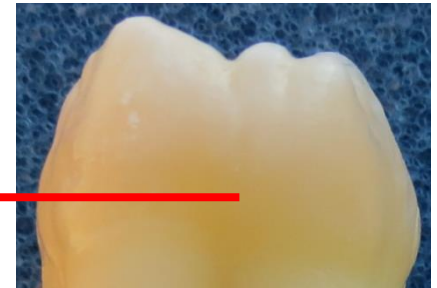

Visually sound
interproximal
surface of Tooth #2

Materials and Methods – The Canary System

Scans from Buccal and Lingual Surfaces

Four areas were scanned at the contact point:


- (1) Distal-buccal (**DB**)
- (2) Distal-lingual (**DL**)
- (3) Mesial-lingual (**ML**)
- (4) Mesial-buccal (**MB**)



Materials and Methods – The Canary System Scans from Occlusal Surfaces

Two areas were scanned at the contact point:

- (1) Distal ridge
- (2) Mesial ridge

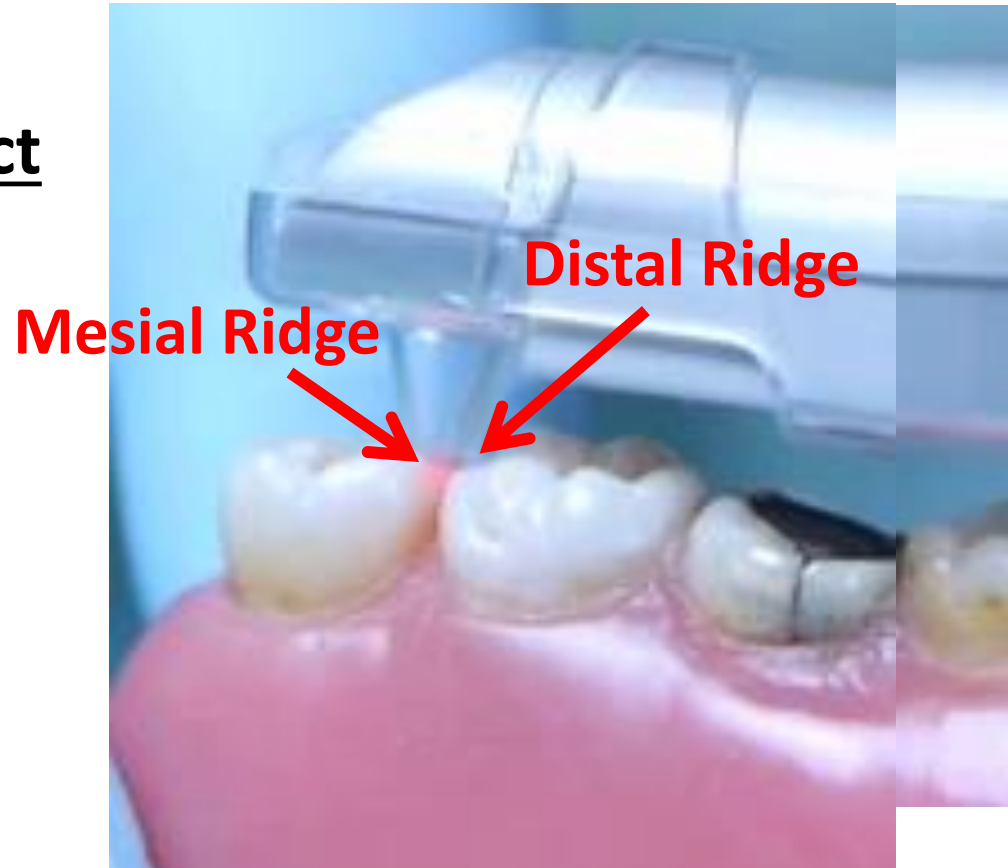

Materials and Methods – ICDAS II and Radiographic Examination

- Two blinded dental clinicians independently scored the interproximal surface of each tooth using:
 1. ICDAS II
 - I. Buccal and lingual
 - II. Occlusal
 2. Radiographs




Materials and Methods - Validation


- Sensitivity and specificity values were calculated according to:
 - i. Direct visual inspection
 - ii. Polarized Light Microscopy (PLM) was performed blinded at the Department of Comprehensive Dentistry, University of Texas Health Science Center at San Antonio as validation



Sensitivity and Specificity – Direct Visual Inspection

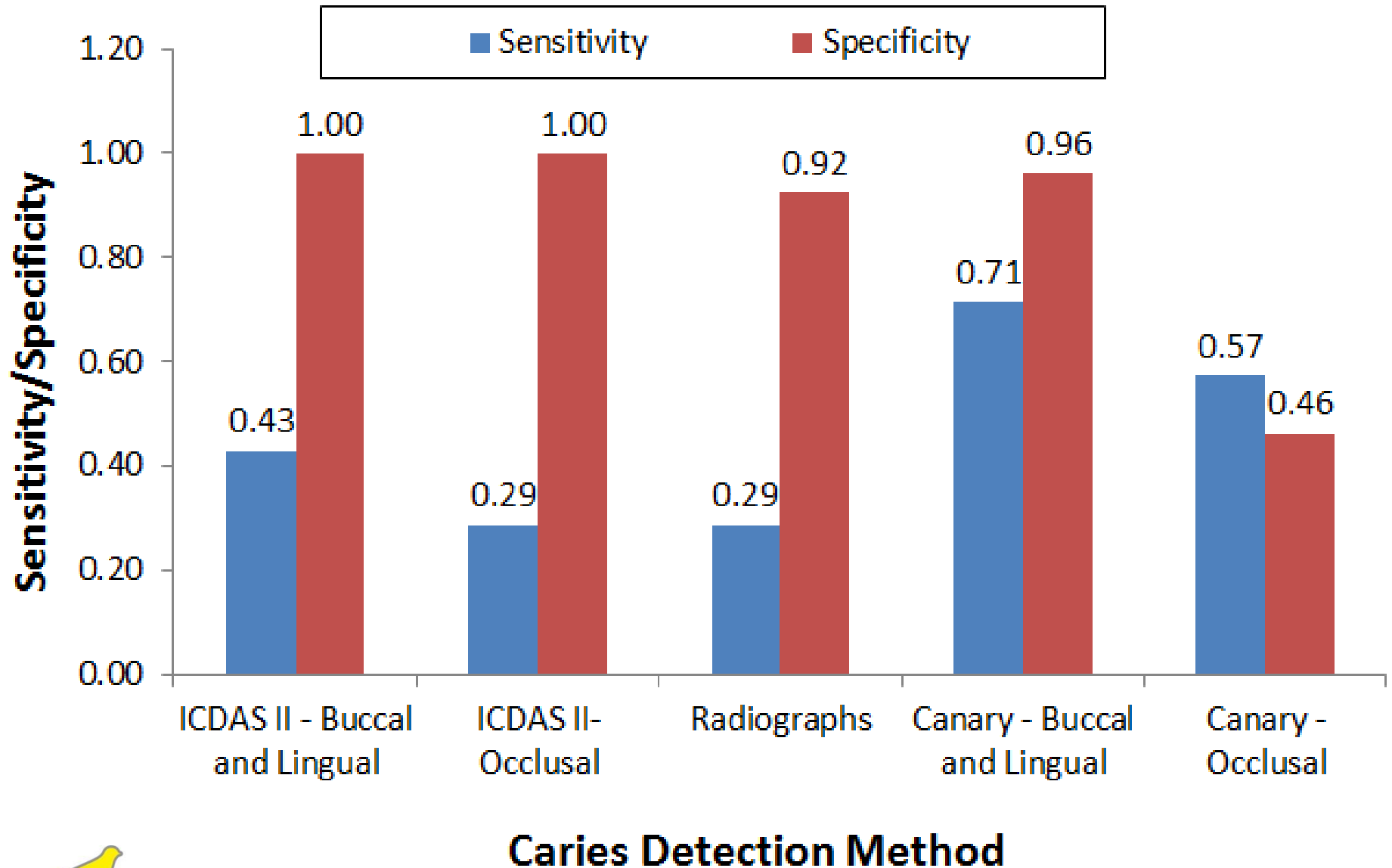

Sensitivity and Specificity – PLM Findings

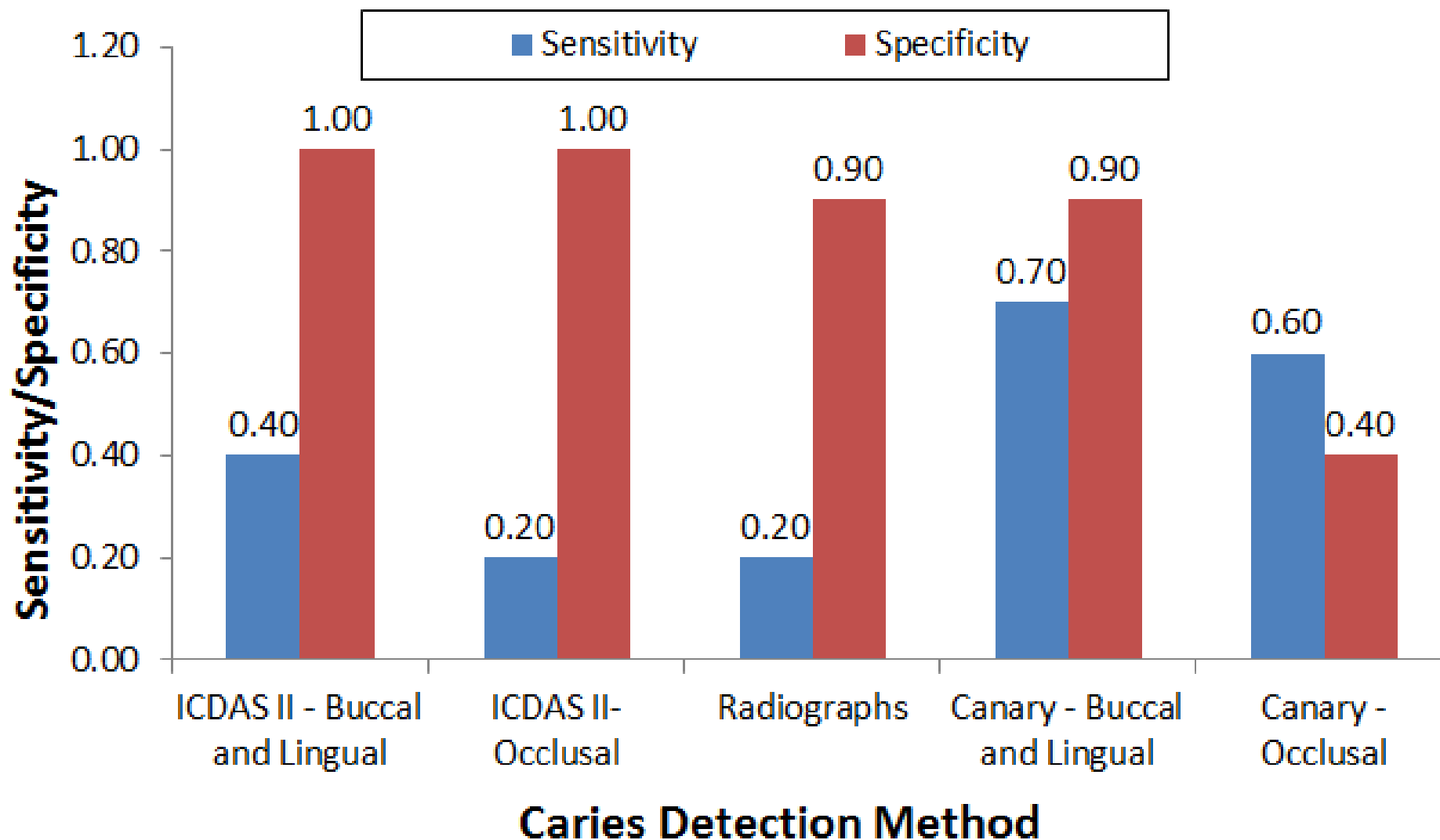


Radiographic Examination

- In this study:
 - High **specificity** (correctly identify absence of caries)
 - Low **sensitivity** (correctly identify presence of caries)
- Similar results reported by previous studies:
 - Pretty IA. J Can Dent Assoc. 2004; 70(6):388–94.
 - Maia AMA e al. Dentomaxillofac Radiol. 2011 Oct;40(7):429-33.
 - Senel B et al. Dentomaxillofac Radiol. 2010 Dec;39(8):501-11.
 - Abesi F et al. Iran J Radiol. 2012 Mar;9(1):17-21.
 - Abreu M et al. Oral Surg Oral Med Oral Pathol Oral Radiol Endod. 2001;91(3):381-5.
 - da Silva Neto JM. Braz Dent J. 2008;19(2):97-102.
 - Bader JD et al. J Dent Educ. 2001 Oct;65(10):960-8.
 - Dove SB. J Dent Educ. 2001;65(10):985-990.
 - Yang J et al. Dent Clin North Am. 2005; 49: 739–752.
 - Pereira AC et al. Caries Res. 2001;35:83-89.
 - Wenzel A. Oral radiology. Principles and interpretation. 5th ed. St. Louis: Mosby; 2004. p. 297–313.



Sensitivity and Specificity – PLM Findings

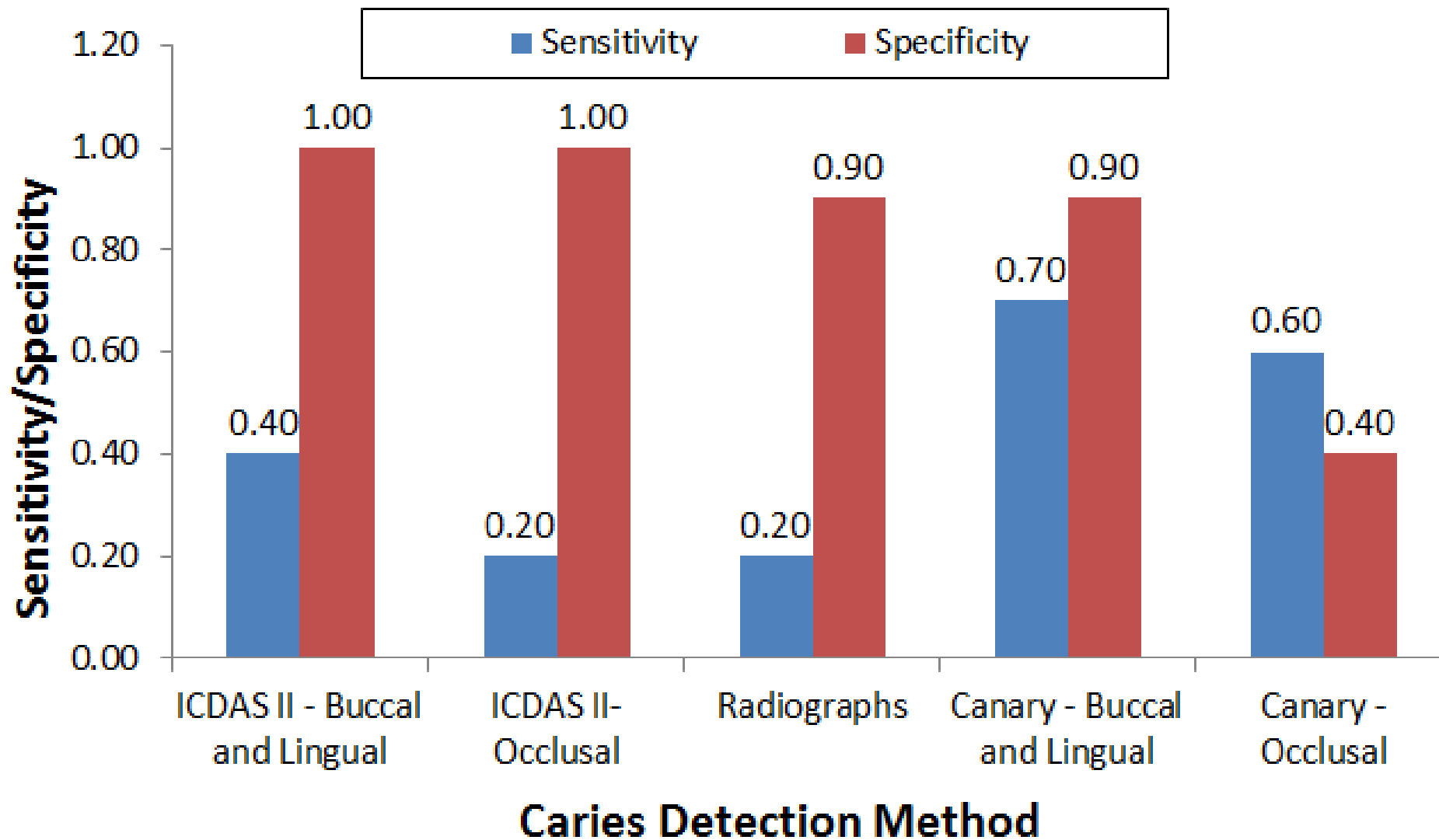

The Canary System

Buccal and Lingual Scans



Conclusions

- This pilot study demonstrated the potential of The Canary System to detect interproximal caries with high sensitivity and specificity when scanning from the buccal and lingual surfaces.
- ICDAS II and radiographic examination resulted in high specificity but poor sensitivity in detecting interproximal caries as previously reported.
- Further *in vitro* studies with larger sample size should be designed to investigate the accuracy and reliability of The Canary System for non-cavitated interproximal caries detection.
- *In vivo* study evaluating the ability of ICDAS II, radiographs, and The Canary System for interproximal caries detection has recently been completed.



The Canary System

by Quantum Dental
Technologies

Thank You

For more information contact:

Bonny Wong

Research Associate, Quantum Dental Technologies

t. (647) 455-2088

e. bonny@thecanarysystem.com